

SIRONIX

Healthium 

REGAIN POWER
WITH SHOULDER IMPLANT
RANGE*



MAY THE SPORT OF LIFE NEVER STOP®


SHOULDER IMPLANTS

BANKART & ROTATOR CUFF REPAIR SOLUTIONS

STATIV[®]

KNOTTED UHMWPE SUTURE ANCHOR

- A soft anchor with patented design provides benefits of bone preservation while giving fixation benefits of traditional anchors.
- Constructed of strong Ultra High Molecular Weight Polyethylene (UHMWPE) suture material.
- This soft anchor is loaded with sutures(s)/tape(s) pre-assembled and loaded on a minimalistic insertion device with bead to give assured centralized deployment of the anchor into pilot hole.
- The deployed unique tri-pod bunched up pattern doubles up compared to the original size of the anchor hence gives much higher pull out strength.

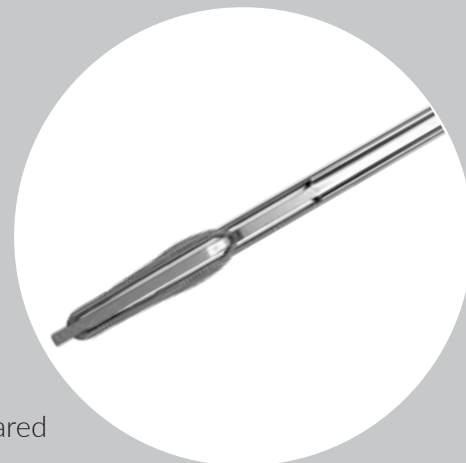

Available Options

Size	1.5mm	1.8mm	2.5mm	2.9mm*
Suture Loading	Single Suture	Double Suture	Double/Triple Suture	Double/Triple Suture
Tape Loading	Single Tape	Double Tape	Double Tape	Double Tape

*Self Tap

VIPLOK[®]

KNOTLESS ANCHOR (PLDLA+ β -TCP / PEEK)

- A versatile knotless anchor comes preloaded on an inserter assembly.
- Intended for fixation of suture to bone/tissue fixation.
- Capable of taking multiple (4-6) sutures/tapes from the tip eyelets*.
- Anchors available in wedge & screw design.



Available Options

Size	2.4mm	4.8mm	5.5mm
Tread Type	Wedge	Screw & Wedge	Screw & Wedge
Anchor Material Type	β TCP & PEEK	β TCP & PEEK	β TCP & PEEK
Tip Material Type	PEEK & Titanium	PEEK & Titanium	PEEK & Titanium


SHOULDER IMPLANTS

BANKART & ROTATOR CUFF REPAIR SOLUTIONS

CEPTRE®

KNOTTED UHMWPE SUTURE ANCHOR (PLDLA+ β -TCP / PEEK / TITANIUM)

- A synergistic combination of unique design, variety of material choices and high strength of UHMWPE suture.
- The unique pattern of threads allows easy insertion and provides maximum fixation strength.
- The anchors are double or triple loaded with suture(s)/tape(s) to provide the best possible combination for superior repair strength.
- The titanium anchor with fully threaded design suitable especially for weak/osteoporotic bone.

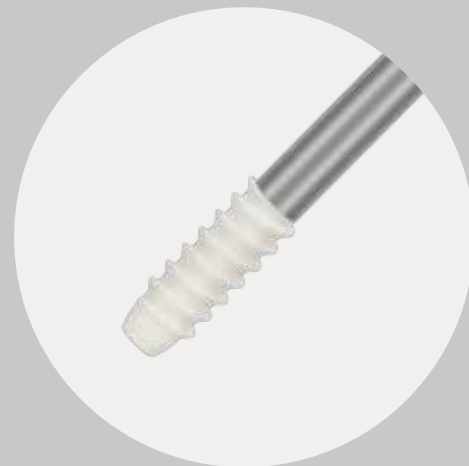

Available Options

Size	4.8mm	5.5mm	6.5mm
Loading Type	Double Suture & Double Tape	Double Suture, Double Tape & Triple Suture	Double Suture, Double Tape & Triple Suture
Material Type	β TCP, PEEK & Titanium	β TCP, PEEK & Titanium	Titanium
	Double Suture with Needle	Suture with Needle	

SPYKE™

KNOTTED ANCHOR

- A shoulder anchor range to offer easy insertion and optimum fixation for instability repair.
- Wedge design allows the cancellous bone to get proper fixation once inserted fully, preventing pullback of anchors.
- The anchors are single or double loaded with suture(s) to provide the best possible combination for superior repair strength.
- The titanium anchor with fully threaded design suitable especially for weak/osteoporotic bone.



Available Options

Size	2.5mm	2.8mm
Loading Type	Single Suture	Double Suture
Material Type	β TCP, PEEK & Titanium	β TCP, PEEK & Titanium

SHOULDER INSTRUMENTS

S36-0001-T	Tissue grasper with Ratchet Straight	S36-2000-T	R-Suture Passer 17mm Working Length
S36-0003-T	Crochet Hook Straight	S39-0011-T	Cannula Introducer (8mm)
S36-0004-T	Switching Stick	S39-0012-T	Cannula Introducer (6mm)
S36-0005-T	Calibrated Probe	S39-2001-T	SP Needle With Lateral Opening
S36-0006-T	Knot Pusher Full Loop	S39-3000-T	SUTURE PASSER-SHOULDER
S36-0007-T	Sliding Suture Cutter Straight	S39-3001-T	Clean Pass Suture Passing Handpiece
S36-0008-T	Clever Hook Right	S39-3002-T	Clean Pass Suture Passing Needle-Medium
S36-0009-T	Clever Hook Left	S39-3003-T	Clean Pass Suture Passing Needle 45 degree LT
S36-0010-T	Sixer Left	S39-3004-T	Clean Pass Suture Passing Needle 45 degree RT
S36-0011-T	Sixer Right	S39-3005-T	Clean Pass Suture Passing Needle 90 degree UP
S36-0012-T	Penetrating Grasper Straight	S39-3006-T	Clean Pass SI Roller replacement Kit
S36-0013-T	Rasp up bend 20 degree	S39-3007-T	Clean Pass Suture Shuttle
S36-0014-T	Tissue Liberator Blade Up	S36-3000-T	Crescent Disposable Suture Passer
S36-0018-T	Tape Cutter	S36-3001-T	Crescent Disposable Suture Passer Needle
S36-0019-T	Suture Manipulator	S36-3002-T	Disposable Lasso Straight
S36-0020-T	Mallet	S36-3003-T	Disposable Lasso 25 Degree Left
S36-0021-T	Shoulder Instrument Box-Large	S36-3004-T	Disposable Lasso 25 Degree Right
S36-0022-T	Shoulder Instrument Set Box-Premium	S36-3005-T	Disposable Lasso 45 Degree Left
S36-0023-T	Combo Grasper	S36-3006-T	Disposable Lasso 45 Degree Right
S36-0024-T	Rasp down bend 20 degree	S36-3007-T	Disposable Lasso 90 Degree Left
S36-0025-T	TissueLiberator Blade Down	S36-3008-T	Disposable Lasso 90 Degree Right
S36-1006-T	Tape knot pusher		

Reference: *Data on file

Disclaimer:

Responsibility for administering products mentioned and interpretation of clinical data provided lies with the registered medical practitioner.

For more information about products, their application & precautions: please visit: www.healthiummedtech.com

Healthium Medtech Limited

RMZ North Star, Cowrks 12th Floor, Adjacent to RMZ Galleria Mall, Yelahanka, Bangalore - 560 064, India
 +(91)-(80)-4186 8198 | healthiummedtech.com