

**SIRONIX**  
Arthroscopy Solutions

**REGAIN AGILITY**  
WITH KNEE IMPLANT  
RANGE\*




**MAY THE SPORT OF LIFE NEVER STOP**

# KNEE IMPLANTS

## FEMORAL FIXATION

### PROLOOP™ (ULTRA, MINI, XL) ADJUSTABLE LOOP BUTTON

- A unique state of the art adjustable loop button fixation device.
- Loop made of non-absorbable, Ultra-High Molecular Weight Polyethylene (UHMWPE) braided durable suture, designed to form an adjustable loop giving surgeons the required adjustability of loop length.
- Proloop has a secured double cross locking mechanism for assured fixation.
- The construct has different colour suture on each side of button for pulling and flipping.

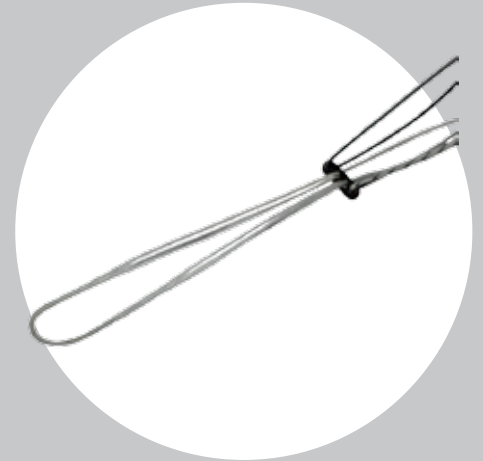

PROLOOP™ ULTRA	PROLOOP™ MINI	PROLOOP™ XL
14mm button with 60mm / 90mm adjustable loop	8mm button with 100mm adjustable loop	20mm button with 90mm adjustable loop

### INFILOOP® FIXED LOOP BUTTON

- A fixed loop button that comes with a fixed length loop and titanium button to provide strong, easy fixation at cortical femur bone.
- The state of the art loop is made of Ultra-High Molecular Weight Polyethylene (UHMWPE) that is durable and provides predictable anchorage of the graft.
- Infiloop is pre-stretched to ensure there is minimal elongation post implantation.



INFILOOP® 12 mm	INFILOOP® 14 mm
Loop length: 10mm, 15mm, 20mm, 25mm, 30mm, 35mm	Loop length: 10mm, 15mm, 20mm, 25mm, 30mm, 35mm


# KNEE IMPLANTS

## TIBIAL FIXATION

### HELYSIS® βTCP/PLLA

- An interference screw provides excellent fixation due to its state of the art patented variable pitch.
- Helysis blunt thread profile & tapered tip design, provides easy insertion and better engagement with graft.
- The design allows wider cannulation via 2mm guide wire and hence avoids tunnel deviation.

βTCP	PLLA
Length 7mm - 11mm Width 20mm-35mm	Length 7mm - 11mm Width 20mm-35mm

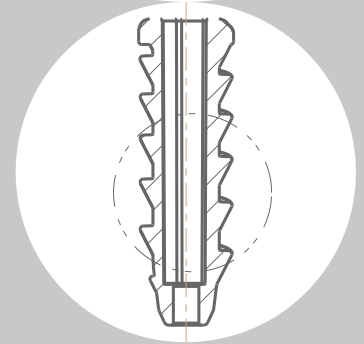

### SURESTITCH™ ALL INSIDE MENISCUS REPAIR IMPLANT

- 17 Gauge stiff needle for effortless piercing\*.
- A built-in adjustable depth control sheath prevents over insertion.
- Safety Knob prevents misfiring and implant slippage.
- Audible & visual confirmation ensures sequential implant selection.
- Ergonomic design enables single handed use.



## OTHER KNEE IMPLANTS

S24-1437-T	Suture Disc D hole Titanium
S24-1438-T	Suture Disc Round hole Titanium
S24-1439-T	Bone Ligament staple (small)
S24-1440-T	Bone Ligament staple (medium)
S24-1441-T	Bone Ligament staple (Large)
S24-1442-T	Spiked washer 4.0mm
S24-1443-T	Cannulated Cancellous Screw 4mmx30mm Short
S24-1444-T	Cannulated Cancellous Screw 4mmx40mm Short
S24-1445-T	Cannulated Cancellous Screw 6.5mmx35mm Short
S24-1446-T	Cannulated Cancellous Screw 4mmx35mm Short
S24-1447-T	Suture Post 30mm

## KNEE INSTRUMENTS

S25-0002-A	Tendon stripper - Open	S25-240D-A	2.4 mm drill bit
S25-0003-A	Tendon stripper - Closed	S25-325D-A	Spade Tip Drill -3.25mm
S25-0004-A	Knee probe	S25-250D-A	Spade Tip Drill - 4.5mm
S25-0005-A	PCL Elevator	S25-1200/1500-A	Guidewire 1.2/1.5 mm
S25-0006-A	Depth Guage	S25-5000/6000/7000/ 8000/9000/1000/1100-A	Cannulated Tibial 5.0/6.0/7.0/8.0/ 9.0/10.0/11.0 mm reamer
S25-0016-A	Graft gauge/ Sizing block	S25-45FD-A/500D/600D/ 700D/750D/800D/850D/ 900D/950D/100D/ 105D/1100	Cannulated Femoral 4.5/5.0/6.0/ 7.0/7.5/8.0/8.5/9.0/9.5/10.0/ 10.5/11.0 mm reamer
S25-0017-A	Graft master	S25-BSNL-A	Helysis Bio Screw driver large
S25-0020/21/22-A	Micro Fracture Awl - 30/45/90 Degree	S25-BSNS-A	Helysis Bio Screw driver Small
S25-0023-A	Staple Impactor	S25-MSDL-A	Screw Driver Fixed-Metal
S25-0024-A	Taper Peg - 7-8mm	S25-0044-A	Cannulated T-handle
S25-0025/26-A	Taper Peg - 8-9/9-10mm	S25-UTDC-A	Sironix Screw Starter
S25-0010/0031/0032-A	Rasp - Straight/Convex/Half Round	S26-0015-T	Shoveler NR Straight Punch
S25-0036-A	Curette	S26-0016-T	Diver Straight Cut Punch
S25-0046-A	Graft Master with Tensioner	S26-0017-T	Diver 15 degree UP
S25-1000-J	ACL Jig Complete Set	S26-0018-T	Grasper Diver With Lock Straight
S25-1000-JT	ACL Jig Complete Set(T)	S26-0019-T	Diver 20 degree Right
S25-1005/1006/1007-A	Femoral Aimer - 5/6/7 mm	S26-0020-T	Diver 20 degree Left
S25-0029-A	Passing Pin		

Reference: \*Data on file

**Disclaimer:**

Responsibility for administering products mentioned and interpretation of clinical data provided lies with the registered medical practitioner.  
For more information about products, their application & precautions: please visit: [www.healthiummedtech.com](http://www.healthiummedtech.com)

**Healthium Medtech Pvt Ltd.**

RMZ North Star, Cowrks 12<sup>th</sup> Floor, Adjacent to RMZ Galleria Mall, Yelahanka, Bangalore - 560 064, India  
+(91)-(80)-4186 8198 | [healthiummedtech.com](http://healthiummedtech.com)